



DEPARTMENT OF
ECOLOGY
State of Washington

Flame Retardants in General Consumer and Children's Products

June 2014

Publication No. 14-04-021

Publication and Contact Information

This report is available on the Department of Ecology's website at <https://fortress.wa.gov/ecy/publications/SummaryPages/1404021.html>

For more information contact:

Hazardous Waste and Toxics Reduction Program

P.O. Box 47600

Olympia, WA 98504-7600

Phone: 360-407-6700

Washington State Department of Ecology - www.ecy.wa.gov

- Headquarters, Olympia 360-407-6000
- Northwest Regional Office, Bellevue 425-649-7000
- Southwest Regional Office, Olympia 360-407-6300
- Central Regional Office, Yakima 509-575-2490
- Eastern Regional Office, Spokane 509-329-3400

Any use of product or firm names in this publication is for descriptive purposes only and does not imply endorsement by the author or the Department of Ecology.

If you need this document in a format for the visually impaired, call the Hazardous Waste and Toxics Reduction Program at 360-407-6700. Persons with hearing loss can call 711 for Washington Relay Service. Persons with a speech disability can call 877-833-6341.

Flame Retardants in General Consumer and Children's Products

by

Saskia van Bergen and Alex Stone

Hazardous Waste and Toxics Reduction Program
Washington State Department of Ecology
Olympia, Washington 98504-7600

Acknowledgements

The authors of this report thank the following Washington State Department of Ecology staff for their contributions to this study:

- Staff at the Department of Ecology's Manchester Laboratory for laboratory analysis, data quality reviews and contract management: Joel Bird, Karin Feddersen, John Weakland, and others.
- ALS Environmental staff for performing the cryomilling of samples and for laboratory analysis.
- The following Department of Ecology headquarters staff:
 - Kelsey Dunne for assistance with X-ray fluorescence (XRF) screening and sample preparation.
 - Josh Grice, Ken Zarker, Carol Kraege and Tom Gries for reviewing the draft report.
 - Cathy Bouge and Linda Riedel for final report editing and publishing.

Table of Contents

Abstract	1
Introduction	3
PBDEs and related compounds	3
Chlorinated phosphate flame retardants	4
Non-halogenated phosphate flame retardants	5
Flame retardants in products	5
Goals and Purpose	6
Methods	6
Sample selection	7
Sample processing and screening	7
Laboratory procedures	9
Data Quality	11
XRF	11
Laboratory	11
PBDEs and deca-BDP-ethane	12
Chlorinated phosphates (from MEL)	12
Non halogenated phosphates (from MEL)	13
HBCD, V6, TBBPA	13
Chlorinated phosphates (from ALS)	13
Non halogenated phosphates (from ALS)	13
Results and Discussion	14
XRF screening	14
Laboratory results from MEL	14
Polybrominated diphenyl ethers and deca-BDP-ethane	15
Chlorinated phosphates	16
Non halogenated phosphates	18
Laboratory results from ALS	19
HBCD, V6, TBBPA	19
Chlorinated phosphates and non halogenated phosphates	19

Summary	22
Recommendations	23
References	24
Acronyms, Abbreviations, and Appendices	27
Appendix 1	28
Appendix 2	35

List of Tables

Table 1.	Number and types of products collected for the 2012-2013 study	7
Table 2.	Analytes of interest	10
Table 3.	Additional analytes of interest	11
Table 4.	Measurement quality objectives for laboratory analyses	12
Table 5.	Bromine XRF screening ranges and matrices sent to MEL greater than 1,000 ppm	14
Table 6.	Summary of matrix and analyte concentration measured in products	14
Table 7.	PBDE and ethane levels measured in products	15
Table 8.	Chlorinated phosphate levels measured in products	17
Table 9.	TPP levels measured in products	18
Table 10.	HBCD, TBBPA and V6 Levels measured in products	19

List of Figures

Figure 1.	General category and distribution of products sent for analysis	8
Figure 2.	General category and distribution of products sent to ALS	9
Figure 3.	Overview of chlorinated phosphate levels in products	16
Figure 4.	Comparison of laboratory levels of TCPP in products	20
Figure 5.	Comparison of laboratory levels of TDCPP in products	20
Figure 6.	Comparison of laboratory levels of TCEP in products	21
Figure 7.	Comparison of laboratory levels of TPP in products	21

Abstract

In 2012-2013, the Washington State Department of Ecology (Ecology) assessed the presence of flame retardant chemicals such as polybrominated diphenyl ethers (PBDEs), in general consumer and children's products.

Ecology collected 169 products from 30 retailers in Washington State between August 2012 and August 2013. Product types included seat cushions, mattresses, upholstered furniture for children, electronics, clothing, and baby carriers. A total of 163 product components from 125 products were sent to a laboratory for analysis.

The majority of samples tested (94%) did not contain PBDEs above a reporting limit of 100 parts per million (ppm), which supports the conclusion that manufacturers have moved away from using PBDEs in products available to Washington State consumers. Alternative flame retardants are still being used. Numerous products were found to contain chlorinated phosphate flame retardants, in particular tris(1,3-dichloro-2-propyl) phosphate (TDCPP) but also tris(2-chloroethyl) phosphate (TCEP) and tris(1-chloro-2-propyl) phosphate (TDPP), and the non-halogenated phosphate triphenyl phosphate (TPP). The majority of these samples were foam and many were children's products.

A subset of samples were tested for TBBPA, HBCD, and a newer flame retardant mixture called Antiblaze® V6 (V6). All three flame retardants were found in some of the samples analyzed.

A number of samples were found to contain a flame retardant identified as a chemical of high concern to children (CHCC) above the reporting limit established in the Children's Safe Product Act (CSPA). Ecology will review the data and determine if the manufacturer is in compliance with the CSPA reporting requirements.

Introduction

PBDEs and related compounds

Polybrominated diphenyl ethers (PBDEs) are a class of persistent, bioaccumulative and toxic (PBT) compounds that historically were used as flame retardants in a wide range of consumer products. In 2001, the total PBDE volume worldwide was estimated at over 67,000 metric tons. There were three different formulations of PBDEs used commercially in products, referred to as penta-BDE, octa-BDE, and deca-BDE. Each was used primarily in different applications. Penta-BDE was used in low density polyurethane foam in products such as furniture. Octa-BDE was used in certain electronic applications such as housings for fax and answering machines, automobile trim, telephone handsets, and kitchen appliance casings. For deca-BDE, the largest use was in electronic enclosures, particularly in computers and televisions (Ecology, 2006). One of the drivers of flame retardants in furniture is the California Standard Technical Bulletin 117, which is a state flammability standard for upholstered furniture.

PBDEs are not chemically bonded to the matrices of the materials. They can escape from their matrix through release to the air and accumulate in dust (Stapleton, 2008).

Studies indicate that PBDEs are ubiquitous throughout the natural environment, in air, soil, and sediments, and are building up in the food chain. PBDEs were detected in migrant Chinook salmon tissue and their stomach contents from four sites in Puget Sound. Other studies have demonstrated thyroid effects on adult fathead minnows and increased risk of disease in juvenile Chinook salmon. There is evidence of bioaccumulation of PBDEs in marine mammals at high concentrations in blubber, including in Puget Sound killer whales. PBDE levels have also been observed to be increasing in orcas resident in the Puget Sound Basin (Ecology, 2006). In conjunction with PBDE concentrations, concerns have been raised about the ‘...*increasing evidence of effects on reproductive health, the immune system, and development in exposed mammals.*’ (Ross, 2006)

Manufacturers of penta-BDE and octa-BDE voluntarily agreed to stop production at the end of 2004 (Ecology, 2006). In 2009, the U.S. Environmental Protection Agency (EPA) announced another voluntary agreement to ‘...*end production, importation and sales of decaBDE for most uses in the United States by December 31, 2012*’ (EPA, 2010). Companies have found alternatives for most PBDE uses.

In addition to the voluntary efforts undertaken by EPA, several states have banned the use of PBDEs in products sold within their state. In 2007, the Washington State Legislature passed legislation ([Chapter 70.76 RCW](#)) banning the use of the penta- and octa-BDE mixtures, banning the use of deca-BDE in one application (mattresses) and providing the possibility of further bans on deca-BDE use if a safer alternative could be identified (Washington, 2007). Ecology and the Washington State Department of Health (DOH) conducted an alternatives assessment to determine if safer alternatives to deca-BDE existed for two additional applications: electronic enclosures (primarily computers and televisions) and residential upholstered furniture. As a result of this assessment, a safer alternative, resorcinol diphenyl phosphate (RDP), was identified in 2009 for electronic enclosures. For upholstered residential furniture, the recommendation was to redesign the product to maintain fire safety without the use of chemical flame

retardants. As of January 1, 2011, the deca-BDE ban was extended in Washington State to include these two applications (Ecology, 2009).

Deca-BDE was also identified by Ecology as a chemical of high concern to children (CHCC). In 2008, the Washington State Legislature passed the Children's Safe Product Act ([Chapter 70.240 RCW](#)) that requires Ecology to identify chemicals that are both toxic and have a potential exposure pathway to children (Washington, 2008). As directed by this legislation, Ecology published a list of 66 CHCCs. Deca-BDE was one of the 66 (Ecology, 2011b). Deca-BDE was placed on this list because of its carcinogenic potential, the tendency of deca-BDE to degrade to more toxic forms of PBDEs such as those found in the penta- and octa-BDE mixtures and children's potential exposure from its wide use as a flame retardant (DOH, 2010). Any manufacturer that uses deca-BDE in children's products sold in Washington State must report the use to Ecology (Ecology, 2011a).

With the banning of PBDEs, manufacturers have switched to a number of alternatives. Brominated alternatives include 2-ethylhexyl 2,3,4,5-tetrabromobenzoate (TBB), (2-ethylhexyl)tetrabromophthalate (TBPH), the polybrominated diphenyl ethanes (PBDP-ethanes), Hexabromocyclododecane (HBCD) and Tetrabromobisphenol-A (TBBPA). Concerns about their toxicological impacts have been raised given similarity to the PBDE formulations. In addition to being CHCCs, HBCD, and TBBPA are also listed as persistent, bioaccumulative toxins (PBTs) by rules adopted in 2006 ([WAC 173-303](#)).

Chlorinated phosphate flame retardants

In addition to deca-BDE, HBDC and TBBPA mentioned above, a chlorinated phosphate flame retardant, tris(2-chloroethyl) phosphate (TCEP), was also identified as a CHCC. Recently, a related compound, tris(1,3-dichloro-2-propyl) phosphate (TDCPP), was identified by the California Environmental Protection Agency as a chemical known to cause cancer and placed on the Proposition 65 list of toxic chemicals (Cal EPA, 2011). It was added to Ecology's CHCC list in 2013 through agency rulemaking. TDCPP has been commonly used as a flame retardant in foam, both in baby products and couches in the United States (Stapleton 2011, 2012). TCEP has also been associated with a newer flame retardant mixture called Antiblaze® V6 (V6), which reports the presence of TCEP as a 4.5-7.5% impurity (EU 2007).

A third flame retardant compound in this class, tris(1-chloro-2-propyl) phosphate (TCPP) has also generated some concern. The Organization of Economic and Cooperative Development (OECD) has indicated that TCPP is harmful to aquatic organisms (OECD, 2012). Given the similarity of this compound to TCEP and TDCPP, concerns have arisen about its widespread use. Like TCEP, TCPP has been associated with a previously uncharacterized flame retardant "U-OPFR" (Stapleton 2011).

Chlorinated phosphates including TCEP have been found in a number of environmental media, such as effluent from sewage treatment plants, precipitation (rain, snow, and glacial ice), and surface waters, primarily in Europe (EU, 2009).

TCEP and TCPP are slightly toxic to aquatic organisms at all trophic levels and TDCPP is moderately toxic to fish. All three compounds are slightly toxic to terrestrial species and aquatic green algae, and are

non-toxic to sewage bacteria (NICNAS, 2001). Organophosphates including TCEP and TDCP have also been found in fish (perch and carp) in a limited study in Sweden (Marklund, 2005).

Non-halogenated phosphate flame retardants

As mentioned previously, Ecology identified RDP as a safer alternative to deca-BDE, especially in the computer and television industry (Ecology, 2009). Other non-halogenated flame retardants were also reviewed in the same report. One of the flame retardants evaluated as an alternative, triphenyl phosphate (TPP), was eliminated due to its high environmental toxicity (Ecology, 2009). Data is lacking on whether these alternative chemicals are being widely used and if concerns about TPP warrant further evaluation.

TPP has been found in:

- Thirty-three dust samples tested in a study in Belgium (Van den Eede, 2011).
- Ninety-eight percent of 50 house dust samples tested in a U.S. study (Stapleton, 2009).
- Fish in a Swedish study of lakes and coastal regions (Sundkvist, 2010).
- Almost half of the fish sampled from Manila Bay in the Philippines (Kim, 2011).
- Reclaimed wastewater (Loraine, 2006).
- The blubber of six bottlenose dolphins from the Gulf of Mexico (Kuehl, 1995).

Flame retardants in products

Recent studies have shown that many of the above flame retardants are found in products and with the phase-out of PBDEs, a more varied group of flame retardants is being used. Stapleton et al. studied polyurethane foam used in baby products and identified flame retardants in a wide range of products including car seats, changing table pads, sleep positioners, portable mattresses, nursing pillows, baby carriers, high chairs, and infant bath mat/slings (Stapleton, 2011). The most common flame retardant used was TDCPP, followed by the mixture Firemaster® 550 (FM 550), which contains multiple chemicals including TBB, TBPH, and TPP.

Chlorinated organophosphate flame retardants were the most common class found including TCEP, TCPP, V6, and “U-OPFR” in addition to TDCPP. Penta-BDE was found in a few of the products with TPP but the authors cautioned these were likely older products manufactured before specific bans were implemented. TPP was also found in one non-halogenated mixture sold commercially as AC073. Excluding demonstrated increase in the number of non-halogenated mixtures, Stapleton et al. found similar results in a study testing couch samples. Penta-BDE was successfully phased-out in the products tested after 2005 (Stapleton, 2012).

The concentration of flame retardants in house dust is widely used to monitor the presence of flame retardants in consumer products. Because many household products contain flame retardants that can be released and accumulated in house dust, dust has become a common sampling media. This connection was investigated and shown to have a high correlation between the amount and type of consumer products and the amount of PBDEs found in house dust (Allen, 2008). Alternative and new flame retardants like HBCD, TBB, TBPH, TCEP, TDCPP, TCPP, V6, and decabromodiphenyl ethane (deca-

BDP-ethane) have been detected in house dust indicating that increasing numbers and different types of flame retardant are now used (Stapleton, 2008, 2009, Fang 2013).

During creation of the CHCC list, Ecology reviewed several sources of data on chemicals in products. Studies conducted by the Danish Environmental Protection Agency and similar bodies of the Dutch Government were particularly useful. In both instances, the organizations purchased products and evaluated them for a wide range of chemicals. This information was used to identify potential sources of exposure to children. Many of the flame retardants in these studies were found in general consumer and children's products. For example, TCEP was found in electronics, perfume in toys and children's articles, and plastic toys. TPP was found in electronics, baby products, and plastic toys. Penta-, octa- and deca-BDE were all found in electronics while deca-BDE was also found in indoor air sampling from various consumer products (Stone, 2010).

Goals and Purpose

The objectives of this project were to evaluate presence of flame retardant chemicals such as PBDEs, in general consumer and children's products and to:

- Determine compliance with Washington State's ban on the PBDE class of flame retardants (RCW 70.76)
- Evaluate the level of substitute flame retardants in various consumer products.
- Verify compliance with Washington's Children's Safe Product Act (CSPA) reporting requirements for flame retardants on the CHCC list (RCW 70.240).

Methods

The design and methods for this study are described in detail in the Quality Assurance Project Plan or QAPP (Ecology, 2012).

Early in 2013 through the CSPA reporting requirements, Ecology was notified that a product sampled for this study contained TBBPA, another flame retardant on the CHCC list, that was not part of the original group of analytes selected for analysis. As part of the effort to determine compliance and assess levels of all pertinent flame retardants, this project was extended to include additional analyses of a subsample of the products. These analyses included three additional flame retardants and the re-analysis of the chlorinated phosphate and triphenyl phosphate flame retardant compounds already analyzed by Ecology's Manchester Environmental Laboratory (MEL). These re-analyses were done to confirm the ability to measure the original analytes. The design and methods for this study are described in a QAPP Addendum (Ecology, 2013).

Sample selection

Products were purchased between August 2012 and August 2013 from local stores in the south Puget Sound area and online retailers. A total of 169 individual products were collected. Products that historically contained flame retardants were targeted for analysis. Table 1 shows the number of samples in each segment identified from the Global Product Classification Standard.

Table 1: Number and types of products collected for the 2012-2013 Study.

Segment	# of products	% of Products collected
Household/Office Furniture/Furnishings	31	18.3
Baby Care	18	10.7
Home Appliances	17	10.1
Clothing	19	11.2
Toys/Games	12	7.1
Electrical Supplies	9	5.3
Building Products	9	5.3
Communications	8	4.7
Computing	8	4.7
Camping	7	4.1
Storage/Haulage Containers	5	3.0
Audio Visual/Photography	5	3.0
Beauty/Personal Care/Hygiene	4	2.4
Cleaning/Hygiene Products	3	1.8
Automotive	3	1.8
Fuels	2	1.2
Home/Business Safety/Security/Surveillance	2	1.2
Lawn/Garden Supplies	2	1.2
Tobacco/Smoking Accessories	2	1.2
Household/Office Furniture	1	0.6
Stationery/Office Machinery/Occasion Supplies	1	0.6
Pet Care	1	0.6
Total Products	169	

Sample processing and screening

Products were processed at Ecology headquarters where each was photographed and assigned a unique identification number. Product information such as brand, country of manufacture, manufacturer, distributor, or importer was recorded.

Products were deconstructed into individual components and screened using X-ray fluorescence (XRF). XRF is a quick, non-destructive method for obtaining the elemental composition of products. Components were defined as individual pieces or parts of the product containing different colors, functions, and/or material. A total of 385 different components from 169 products were screened for bromine with a scan

time of 120 seconds using a Niton XL3t XRF analyzer. Chlorine was not effectively screened using XRF due to instrument settings. All samples were measured in Niton's TestAll® or Plastics Non PVC mode. Details of the sample preparation and XRF screening are described in the [Sample Processing Standard Operating Procedure](#) (publication 14-04-013). Chain-of-custody was maintained and documented throughout the process.

A total of 163 component samples from 125 products were sent to MEL for analysis. Fabric, stuffing, and foam samples were manually cut into small pieces. Plastic samples and some foam samples were frozen and milled in a cryogenic chamber. Plastic samples that contained metal, for example motherboards, were milled and metal removed. The components, type of products, and material identification (foam, plastic, etc.) were identified in the chain of custody (COC) sent to the laboratory ([Appendix 1](#)). General categories of the products sent to MEL and distributions are shown in [Figure 1](#). Forty percent of the products were targeted towards children including baby activity gyms, car seats, changing mats, pajamas, children's chairs, and stuffed animals. Seventy-two (44.2 %) of these samples were plastic, 42 (25.8 %) were foam, 33 (20.2 %) were fabric, 13 were stuffing (8.0%), one (0.6%) was fiberglass, one (0.6%) was Styrofoam, and one (0.6 %) was a solution.

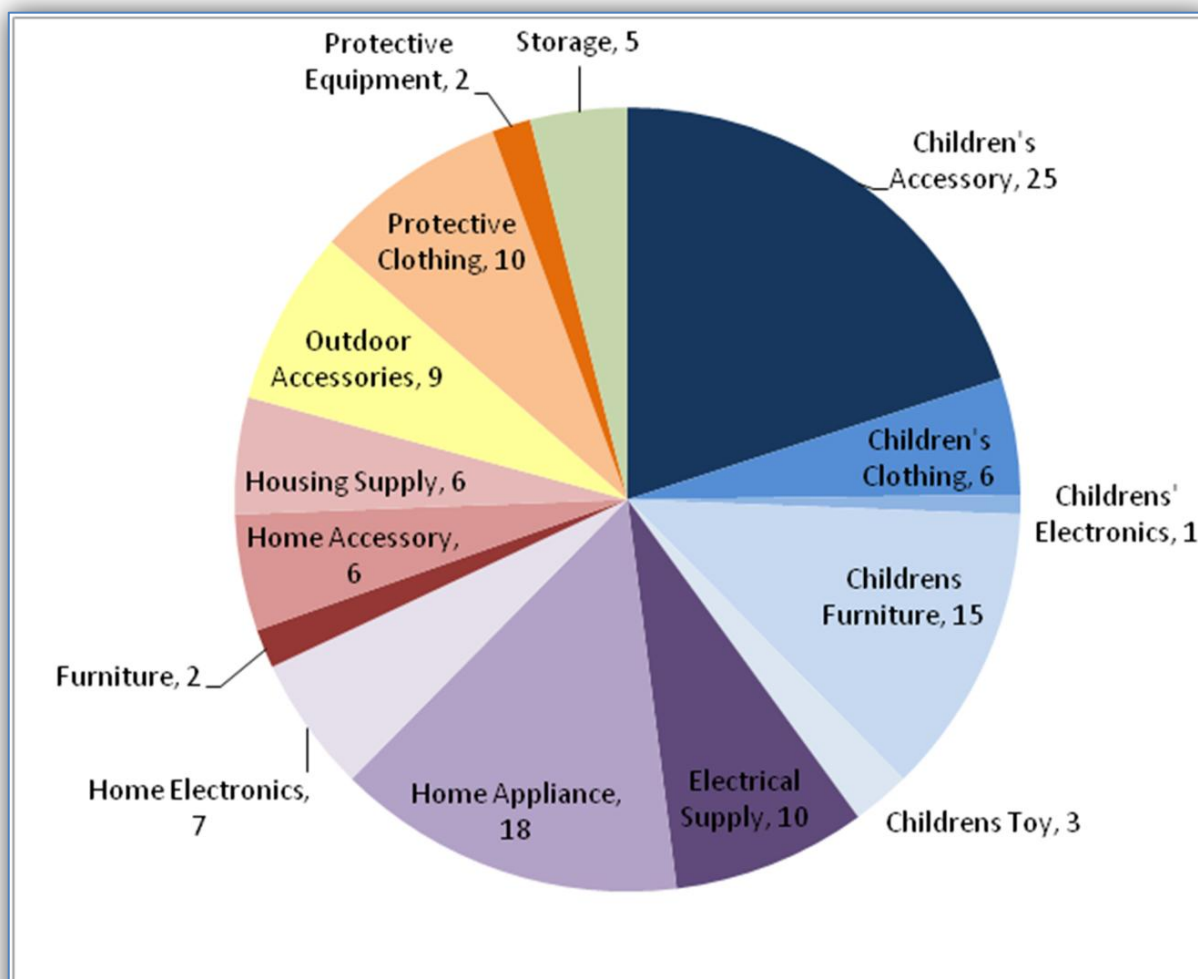


Figure 1: General category and distribution of products sent for analysis.

From these 163 components, a subsample of 67 components from 61 products was sent to a second laboratory (ALS Environmental or ALS) to be analyzed for three additional analytes. The chlorinated and non-halogenated phosphates were re-analyzed to further validate MEL's analytical results. The bolded samples listed in [Appendix 1](#) were sent to the second laboratory. The samples were chosen based primarily on bromine XRF screening results and preliminary data from MEL. Secondly, a sample was sent if it had a California Standard Technical Bulletin 117 label (a state flammability standard for upholstered furniture) as the label indicates that a flame retardant is very likely present (Stapleton 2012). The distribution of products is displayed in [Figure 2](#).

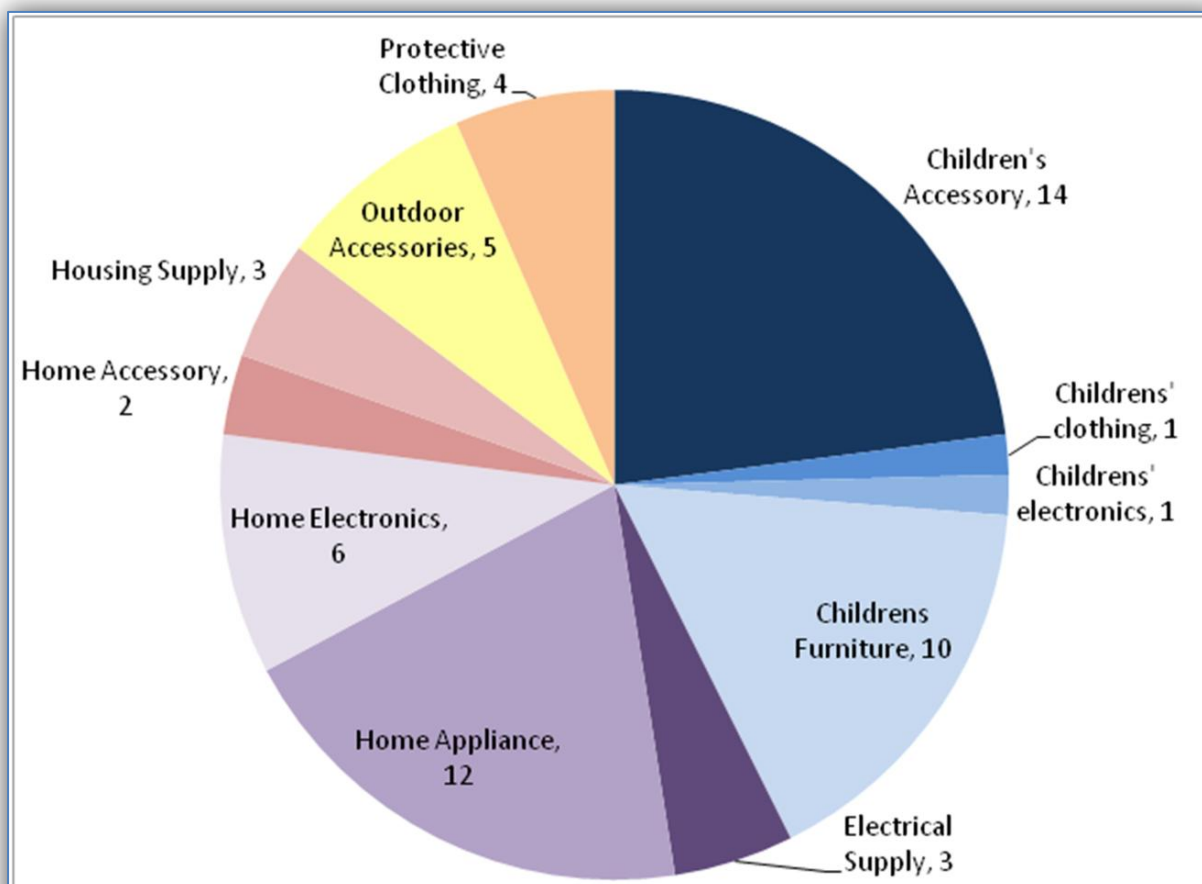


Figure 2: General category and distribution of products sent to ALS for validation.

Laboratory procedures

Cryomilling/Milling

Plastic and some foam samples selected for laboratory analysis were sent to ALS for cryomilling prior to analysis. Samples that could be cryomilled were done following ALS' Standard Operating Procedure for pulverizing samples employing a freezer/mill (EXT-GRIND). Samples were made brittle using cryogenic temperatures and then milled in an enclosed capsule by magnetically shaking the capsule back and forth. After cryomilling, samples consisted of homogenous particles.

Flame Retardant Analysis

MEL extracted the flame retardant compounds listed in Table 2 in all samples following the Environmental Protection Agency's (EPA's) extraction Method 3546 (MEL Standard Operating Procedure (SOP) 730122) and analyzed using gas chromatography-mass spectrometry (GC/MS) following a modified EPA analytical Method 8270D (MEL SOP 730123). RDP (CAS #115-86-6) was included in the QAPP but was not analyzed due to lack of pure standard material.

Table 2. Analytes of interest

Analytes	CAS Number
PBDEs	
PentaBDE formulation:	
PBDE 47	5436-43-1
PBDE 49	243982-82-8
PBDE-66	189084-61-5
PBDE-71	no CAS available
PBDE-99	60348-60-9
PBDE-100	189084-64-8
PBDE-138	182677-30-1
OctaBDE formulation:	
PBDE 153	68631-49-2
PBDE-154	207122-15-4
PBDE-183	207122-16-5
PBDE-184	189084-68-2
PBDE-191	no CAS available
Decabromodiphenyl ether (PBDE 209)	1163-19-5
Decabromodiphenyl ethane (deca-BDP-ethane)	84852-53-9
Chlorinated phosphates	
TCEP	115-96-8
TCPP	13674-84-5
TDCPP	13674-87-8
Non-halogenated phosphates	
TPP	125997-21-9

Approximately 0.20 grams of sample were extracted via microwave-assisted extraction using a 70:30 acetone:hexane solution. After extraction, samples were concentrated, solvent exchanged into isooctane, and diluted. Samples were spiked with surrogates and an internal standard.

A subset of the initial samples was analyzed for the chlorinated and non-halogenated phosphates listed in Table 2 and the additional analytes listed in Table 3. ALS extracted these flame retardant compounds using EPA extraction Method 3540C. The extracts were analyzed using two methods:

1. Gas chromatography-mass spectrometry (GC/MS) following a modified EPA Method 8270D for the analytes in Table 2.
2. Liquid chromatography-mass spectrometry/ mass spectrometry (LC/MS/MS) following a modified EPA Method 1694 for the analytes in Table 3.

Table 3: Additional analytes of interest

Analytes	CAS Number
V6	38051-10-4
HBCD	3194-55-6 25637-99-4
TBBPA	79-94-7

Approximately 0.20 grams of sample was extracted via Soxhlet extraction using a 70:30 acetone:hexane solution. After extraction, samples were concentrated and split. The extract for LC/MS/MS was solvent exchanged into acetonitrile. Samples were spiked with surrogate and internal standard solutions.

Standards were purchased as pure neat material with the exception of V6, which was purchased as a technical grade material and purified following the procedure described in Fang et al. 2013.

Data Quality

XRF

XRF performance was assessed at the beginning of a batch of samples with a system test and the results of a plastic reference sample run at the beginning and end of a batch. The plastic standard did not include the elements bromine and chlorine.

Laboratory

MEL conducted a data quality review of their data packages and the sub-contracted laboratory data packages. Case narratives describing the quality of laboratory data, including instrument calibration, and quality control tests, are available upon request. Quality control tests for each batch analyzed consisted of a method blank, laboratory control samples (LCS), matrix spikes/matrix spike duplicates, and if enough sample was available, a laboratory duplicate. The reporting limit of 100 ppm was selected for the analytes of interest.

Measurement quality objective (MQO) targets outlined in the QAPP and Addendum differed from those reported from the lab as there was minimal history on the methods and the initial target criteria were too restrictive. The lab defined changes to the criteria and the Project Manager agreed. The actual acceptance criteria are shown in Table 4 and reviewed below.

Table 4: MQOs for laboratory analyses

Analyte	Lab Control Samples (recovery)	Matrix ⁺ Spikes (recovery)	Duplicates ⁺ (RPD) ⁺⁺	Method Blanks (ppm) ⁺⁺⁺	Surrogate Recovery (recovery)
PBDEs	50-150%	50-150%	± 40%	<100 for all except <500 for 209	50-150%
Deca-BDP-ethane	50-150%	50-150%	± 40%	< 100	50-150%
Chlorinated phosphates	50-150%	50-150%	± 40%	< 100	50-150%
Non-halogenated phosphates	50-150%	50-150%	± 40%	< 100	50-150%
V6	60-140%	60-140%	± 50%	< 100	50-150%
HBCD	60-140%	60-140%	± 50%	< 100	50-150%
TBBPA	60-140%	60-140%	± 50%	< 100	50-150%

⁺ Matrix spike duplicates and split duplicates

⁺⁺RPD = Relative Percent Difference

⁺⁺⁺ ppm = parts per million

PBDEs and deca-BDP-ethane

Data for the sample analyses were generally within the MQO targets outlined in Table 4. All method blanks were below reporting limits. MEL encountered solubility problems with deca-BDP-ethane analytical standards; therefore, all deca-BDP-ethane results were qualified as estimates. PBDEs and deca-BDP-ethane were not spiked into the Matrix Spikes due to cost considerations.

Instances where MQOs were not achieved or standard laboratory procedures were outside of acceptance limits included:

- Several calibration checks (continuous) for deca-BDP-ethane were outside of acceptance limits resulting in qualifying the data as estimates.
- Several surrogate recoveries were outside of the MQO limits set resulting in qualifying the brominated data as estimates.
- Laboratory Control Sample (LCS) recoveries were low in one BDE-209 batch resulting in a UI flag.

Additionally, instrument calibration procedures were outside of laboratory-identified targets in several instances. This, along with the MQO violations above, caused some results to be qualified as estimates.

Chlorinated phosphates (from MEL)

Data for the sample analyses were generally within the MQO targets outlined in Table 4. All method blanks were below reporting limits. Instances where MQOs were not achieved or standard laboratory procedures were outside of acceptance limits included:

- Several calibration checks (initial) for TCP and TDCPP were outside of acceptance limits resulting in qualifying the data as estimates.

- Several calibration checks (continuous) for TCPP and TCEP were outside of acceptance limits resulting in qualifying the data as estimates.
- An LCS for TCPP was outside acceptance limits resulting in qualifying the data as estimates.
- Recoveries of one matrix spike were not calculated for TCPP and TDCPP due to high concentrations in the native sample.
- Recovery of one matrix spike was low for TCEP resulting in a UJ flag for the native sample.

Non halogenated phosphates (from MEL)

Sample data were generally within the MQO targets outlined in Table 4. All method blanks were below reporting limits. One method blank tentatively identified TPP below the reporting limit but all samples with detects were greater than five times the amount found in the blank; therefore, no qualifiers were added. Instances where MQOs were not achieved or standard laboratory procedures were outside of acceptance limits included:

- Surrogate recoveries of TPP-d15 were outside of acceptance limits in some samples resulting in qualifying the TPP data as estimated.
- Recoveries of one matrix spike were not calculated for TPP due to high concentrations in the native sample.

HBCD, V6, TBBPA

Sample data were within the MQO targets outlined in Table 4 with the exception of one matrix spike that had a high recovery of HBCD. All method blanks were below reporting limits.

Chlorinated phosphates (from ALS)

Data for the sample analysis were generally within the MQO targets outlined in Table 4. All method blanks were below reporting limits. Instances where MQOs were not achieved or standard laboratory procedures were outside of acceptance limits included:

- One matrix spike recovery was not calculated as the native sample was inhomogeneous.
- TDCPP recoveries for one matrix spike were outside of acceptance limits resulting in the native sample being qualified as an estimate.
- One duplicate sample had TDCPP five times lower than the first sample, which resulted in qualifying both samples as estimates.

Non halogenated phosphates (from ALS)

Sample analysis data were generally within the MQO targets outlined in Table 4 with the exception that surrogate recoveries of TPP-d15 were outside of acceptance limits in two samples resulting in qualifying the TPP data as UJ. All method blanks were below reporting limits.

Results and Discussion

XRF screening

Three hundred, eighty-five different components from 169 products were screened for bromine using an XRF analyzer. Chlorine was not effectively screened using XRF due to instrument settings.

Scanning 385 components using XRF produced:

- 207 samples with non-detects for bromine (<LOD)
- 110 samples with bromine levels less than 1,000 ppm
- 26 samples with bromine levels between 1,000 and 5,000 ppm
- 42 samples with bromine levels above 5,000 ppm (0.5%)

Of the samples greater than 1,000 ppm, 79.4% of them (54 samples) were sent to the laboratory. Most of these were plastic. Table 5 lists the matrix type of the samples sent to MEL above 1,000 ppm bromine.

Table 5: Bromine XRF screening ranges and matrices sent to MEL greater than 1,000 ppm

Matrix	Bromine XRF Results	
	1,000-5,000 ppm	> 5,000 ppm
Foam	7	4
Plastic	7	29
Fabric	4	3

Laboratory results from MEL

One hundred, sixty-nine component samples from 125 products were sent to MEL laboratory for analysis of the flame retardants identified in Table 2. Forty (23.7%) of these components from 36 (28.8 %) products had concentrations of one or more of the analytes (Table 2) above the method reporting limit. Individual concentrations ranged from the reporting limits to 84,000 ppm. A summary of the amount of analyte and matrix are listed in Table 6. Most of the samples with analytes in the percent level (>10,000 ppm) were foam samples. The remaining two were fabric.

Table 6: Summary of matrix and analyte concentration measured in products

Matrix	Sum of analytes			Total samples above RL
	< 1000 ppm	1,000 - 9,999 ppm	> 10,000 ppm (> 1%)	
Foam	1	7	17	25
Plastic	7	3	-	10
Fabric	1	2	2	5

Polybrominated diphenyl ethers and deca-BDP-ethane

Table 7 displays the components where the PBDEs and deca-BDP-ethane concentrations measured above the reporting limits. These analytes were detected in ten components from nine products. PBDEs were found in eight components of seven products and deca-BDP-ethane was found in six components of five products. The majority of the samples were very close to their reporting limits. Only three results, which are bolded in Table 7, were above 1,000 ppm. As mentioned previously, MEL encountered solubility problems with deca-PBD-ethane analytical standards; therefore, all deca-BDP-ethane results were qualified as estimates. One sample was from a children's product (product description is highlighted blue).

Table 7. PBDE and ethane levels measured in products

Sample Name	Analyte	Result (ppm)	Matrix	Product Description
GR000-F01	BDE-209	2,400	Plastic	Plastic Pallet
OM003-F01	BDE-183	190	Plastic	Shredder
	BDE-209	2,600		
	Deca-BDP-ethane	750 J		
OM003-F02	BDE-183	310	Plastic	Shredder
	BDE-209	1,800		
	Deca-BDP-ethane	1,000 J		
LW002-F01	BDE-047	150	Foam	Carpet Padding
LW000-F01	BDE-209	200	Foam	Carpet Padding
	BDE-047	150		
	BDE-099	300		
LW001-F01	BDE-209	140	Foam	Carpet Padding
	BDE-047	110 NJ		
	BDE-099	250		
WM066-F01	Deca-BDP-ethane	250 NJ	Plastic	Battery Charger
WM068-F03	Deca-BDP-ethane	470 J	Plastic	Foot Warming Pad
CT009-F03	Deca-BDP-ethane	350 NJ	Plastic	LED TV
BY002-F08	Deca-BDP-ethane	620 J	Plastic	Child's Tablet

In Washington, penta-BDE and octa-BDE are prohibited in all products. Deca-BDE is prohibited in mattresses, residential upholstered furniture, and "electronic enclosures" for computer monitors and televisions. Ecology has issued [guidance](#) clarifying that it may take enforcement actions if PBDE levels exceed 1,000 ppm by weight in a restricted product. None of the PBDE levels found violated Washington law.

Flame retardants are typically used in the percent level. Levels near the reporting limit could be cross contamination during manufacture, flame retardant impurities, or flame retardants from recycled content. Concentrations lower than percent level could be part of a mixture. While 54 of these samples screened positive for bromine above 1,000 ppm using the XRF, only three samples contained brominated analytes listed in Table 2 above 1,000 ppm. This suggests that there are additional brominated compounds in the samples.

Firemaster® 550 mixture (FM550) is one likely alternative to PBDEs. It contains two brominated compounds, 2-ethylhexyl-2,3,4,5-tetrabromobenzoate (TBB) and bis(2-ethylhexyl)-2,3,4,5-tetrabromophthalate (TBPH), and triphenylphosphate (TPP) in addition to several isopropylated triaryl phosphate isomers. FM550 has been found in children's products and couches (Stapleton 2011, 2012). Two additional alternatives are HBCD and TBBPA. A subset of samples from this study was evaluated for the latter two flame retardants and these results will be discussed in a [later section](#).

Chlorinated phosphates

Chlorinated phosphate flame retardants were detected in 27 components (21 foam, five fabric, and one plastic) from 24 products. Twenty-four of the components from 22 products contained individual flame retardants above 1,000 ppm. Of these 22 products, 14 are children's products. About half of the samples contained flame retardants at the percent level (above 10,000 ppm). Table 8 displays the list of samples where at least one chlorinated phosphate level measured above the reporting limit.

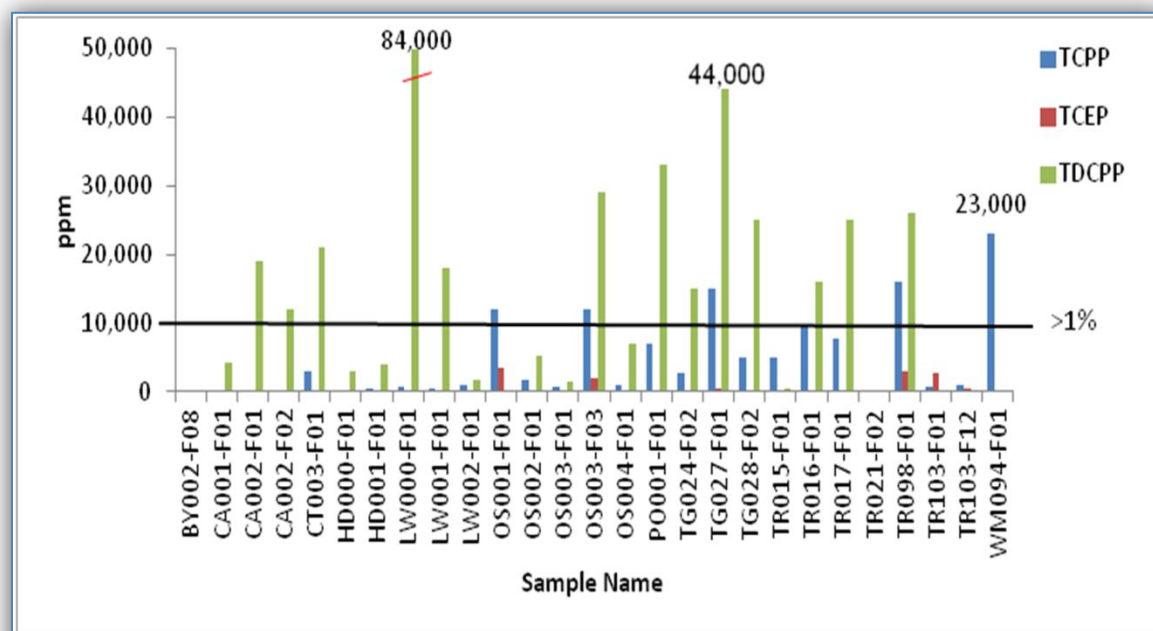


Figure 3: Overview of chlorinated phosphate levels in products.

In agreement with previous studies (Stapleton 2011, 2012), TDCPP was the most common chlorinated phosphate detected in foam. Twenty samples (17 foam and three fabric) contained TDCPP above 1,000 ppm. The fabric components came from tent products, which met the CPAI-84 Flammability Specification for Flame-resistant Materials Used in Camping Tentage, the U.S. tent flame retardant standard promulgated by the Industrial Fabrics Association International. Keller et al (2014) found similar levels in tent textiles.

Four children's products contained TCEP above 1,000 ppm. Seven samples from six children's products were above state Practical Quantitation Limit (PQL) of 100 ppm for reporting TCEP to Ecology as required

by the CSPA. Five samples were foam while two were fabric. These values are bolded in red in Table 8. TCEP can be an impurity from another flame retardant known as V6, also found in foam. A subset of samples from this study was evaluated for V6. The results are discussed in a [later section](#).

Thirteen products, all foam, contained TCPP above 1,000 ppm. Eleven samples were children's products. Many of these also contained a higher concentration of TDCPP. Similar to TCEP, TCPP can be an impurity in another flame retardant found in recent studies currently called 'unidentified chlorinated organophosphate flame retardant' (U-OPFR).

Table 8: Chlorinated phosphate levels measured in products

Sample	Analytes (ppm)			Matrix	Product Description
	TDCPP	TCPP	TCEP		
TG027-F01	44,000 J	15,000 J	550 J	Foam	Booster Seat
OS003-F03	29,000 J	12,000 J	2,100 J	Foam	Child's Chair
TR098-F01	26,000	16,000 J	3,000	Foam	Child's Chair
TG028-F02	25,000 J	4,900 J	< 87	Foam	Changing Pad
TR017-F01	25,000 J	7,700	< 93	Foam	Crib Wedge
TR016-F01	16,000 J	9,700	< 92	Foam	Portable Crib Pad
TG024-F02	15,000	2,800 J	270 J	Foam	Changing Mat
OS004-F01	7,000 J	890	< 98	Foam	Child's Chair Accessory
OS002-F01	5,300	1,800 J	< 97	Foam	Child's Chair
OS003-F01	1,600	700	100 J	Foam	Child's Chair
TR015-F01	370 J	5,000	< 93	Foam	Changing Pad
TR021-F02	270	< 98	< 98	Foam	Pee protector
BY002-F08	250 J	< 92	< 92	Plastic	Child's Tablet
TR103-F12	< 97	870	510 J	Fabric	Baby Carrier
OS001-F01	< 97	12,000	3,400 J	Foam	Child's Chair
TR103-F01	< 96	640 J	2,700	Fabric	Baby Carrier
WM094-F01	< 96	23,000 J	< 96	Foam	Child's Chair
LW000-F01	84,000 J	690	120	Foam	Carpet padding
PO001-F01	33,000 J	7,000 J	< 98	Foam	Seat Cushion
CT003-F01	21,000 J	3,000	< 96	Foam	Chair Pad
CA002-F01	19,000 J	< 100	160	Fabric	Tent
LW001-F01	18,000 J	370	< 97	Foam	Carpet padding
CA002-F02	12,000 J	< 98	130	Fabric	Tent
CA001-F01	4,300	ND	180	Fabric	Tent floor liner
HD001-F01	4,100 J	550	< 94	Foam	Carpet padding
HD000-F01	2,900 J	180	< 96	Foam	Carpet padding
LW002-F01	1,800 J	990	150	Foam	Carpet padding

Blue highlighted results are children's products.

Non halogenated phosphates

Fifteen samples from 14 products contained TPP above the reporting limits (Table 9), where seven were above 1,000 ppm. TPP is used both as a plasticizer and as a flame retardant in halogenated and non-halogenated flame retardant mixtures (Stapleton 2009, 2012). The two non-halogenated mixtures found by Stapleton et al. contained TPP with either tris(4-(tert-butyl)phenyl) phosphate and several butylphenyl isomers or with several methyl- or dimethyl- phenyl phosphate isomers. One halogenated mixture containing both TPP and the halogenated compounds 2-ethylhexyl 2,3,4,5-tetrabromobenzoate (TBB) and (2-ethylhexyl)tetrabromophthalate (TBPH), is known as FM550. FM550 has frequently been found in children's foam products and furniture (Stapleton 2011, 2012).

Of the seven samples with TPP above 1,000 ppm, four had bromine levels above 10,000 ppm (one percent by weight) suggesting these products used halogenated mixtures. Three were children's furniture containing foam.

Table 9: TPP levels measured in products

Sample	Analyte	Conc. (ppm)	Matrix	Product Description	Br > 1% (XRF)
AM009-F01	TPP	29,000 J	Foam	Child's Chair	
TG031-F01	TPP	13,000 J	Foam	Child's Chair	+
AM008-F01	TPP	10,000 J	Foam	Child's Chair	+
PB000-F01	TPP	6,300	Foam	Child's Chair	+
OM005-F08	TPP	21,000	Plastic	LCD Monitor	+
WM066-F01	TPP	13,000	Plastic	Battery Charger	
HD001-F01	TPP	2,700	Foam	Carpet padding	
CT010-F05	TPP	920	Plastic	Blue ray DVD player	
LW000-F01	TPP	700	Foam	Carpet padding	
HD000-F01	TPP	610	Foam	Carpet padding	
LW002-F01	TPP	560	Foam	Carpet padding	
OM003-F01	TPP	440	Plastic	Shredder	+
LW001-F01	TPP	400	Foam	Carpet padding	
OM003-F02	TPP	290	Plastic	Shredder	+
OM000-F02	TPP	210	Plastic	LCD Monitor	

+ samples screened positive for > 1% bromine using XRF

Laboratory results from ALS

A subsample of Sixty-seven component samples from 61 products were sent to a second laboratory (ALS) to be analyzed for three additional analytes listed in Table 3, and the chlorinated and non-halogenated phosphates from Table 2 were re-analyzed to compare the results from the two laboratories.

HBCD, V6, TBBPA

HBCD, V6, and TBBPA were detected in 14 samples from 14 products (Table 10). Nine were above 1,000 ppm.

Table 10: HBCD, TBBPA, and V6 levels measured in products

Sample	Analyte	Conc. (ppm)	Matrix	Product Description
GR003-F01	HBCD	44,300	Plastic	Protective Glove
WM095-F01	HBCD	563	Styrofoam	Bean Bag Chair
WS002-F01	TBBPA	94,300	Plastic	Car charger
WM068-F03	TBBPA	65,100	Plastic	Heated mattress foot warming pad
OM003-F02	TBBPA	35,200	Plastic	Shredder
OM000-F02	TBBPA	12,400	Plastic	Battery Charger
BL005-F01	TBBPA	282	Plastic	Hair Flat Iron
FM043-F01	TBBPA	196	Plastic	Clothing Iron
TR103-F01	V6	42,500	Foam	Baby Carrier
OS001-F01	V6	37,200	Foam	Child's Chair
TR098-F01	V6	16,200	Foam	Child's Chair
CA002-F01	V6	2,390	Fabric	Tent
LW002-F01	V6	201	Foam	Carpet padding
LW000-F01	V6	74.5	Foam	Carpet padding

HBCD was detected in a plastic and a Styrofoam component. One sample was above 1,000 ppm and the other, a child's bean bag chair, was above the CSPA criteria for reporting HBCD in children's products. This potential compliance value is bolded in red (Table 10).

V6 was detected in one fabric and five foam products. Four of these results were above 1,000 ppm. TCEP is an impurity in V6 and has been reported as 4.5-7.5 % of the V6 (EU, 2007). Three of these children's products were above the CSPA criteria for TCEP in children's products, which indicates that V6 was the source of TCEP.

TBBPA was detected in six plastic components in six products. Four of these were above 1,000 ppm. None of these were children's products.

Chlorinated phosphates and non halogenated phosphates

Samples containing flame retardants above 1,000 ppm from ALS or MEL are listed in [Appendix 2](#). Thirty-nine analytes were detected in 23 different products and the majority were foam components. Detection consistence between laboratories was excellent indicating that MEL and ALS were able to measure the analytes reliably. Only a qualitative comparison is possible as no reference standards in the appropriate matrices and only two labs were used.

The results (Figures 4-7) indicate ALS reported higher values than MEL in almost all samples. This would indicate that either ALS is biased high or MEL is biased low or both. Sample LW000-F01's TDCPP result was the only exception. This sample was a foam carpet padding made up of multiple colored foam. It is believed that the sample was not homogeneous and was the main cause for the variable results.

The accuracy of the results from the two labs is unknown as there are no reference standards for the analytes measured in the matrices tested. Some of the main factors that affect accuracy are:

- Extraction efficiency, which can be affected by extraction type, matrix, solvent used and extraction conditions.
- Analytical methods including surrogates and internal standards used, instrumentation and calibration method.

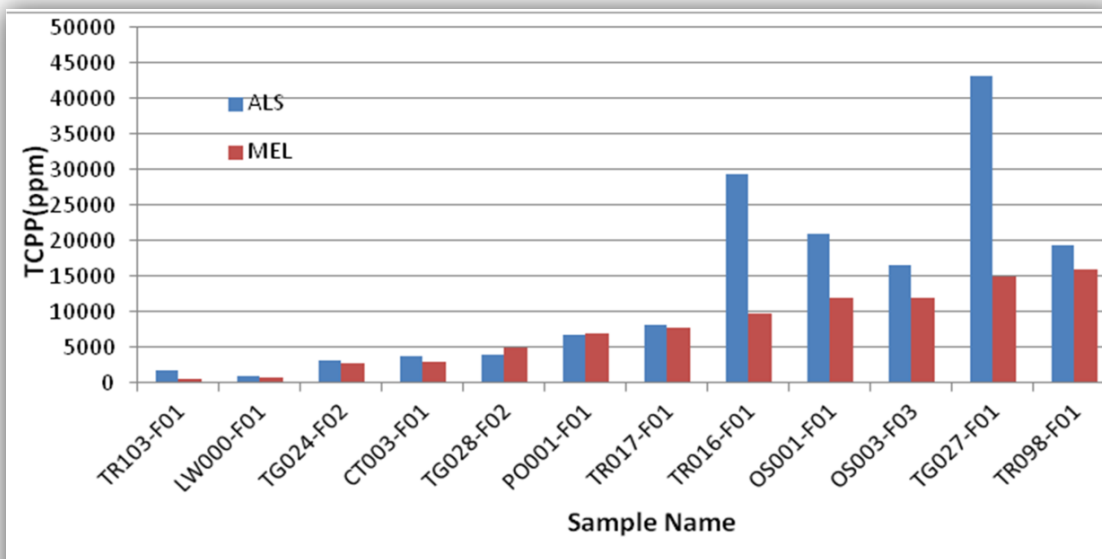


Figure 4: Comparison of laboratory levels of TCPP in products.

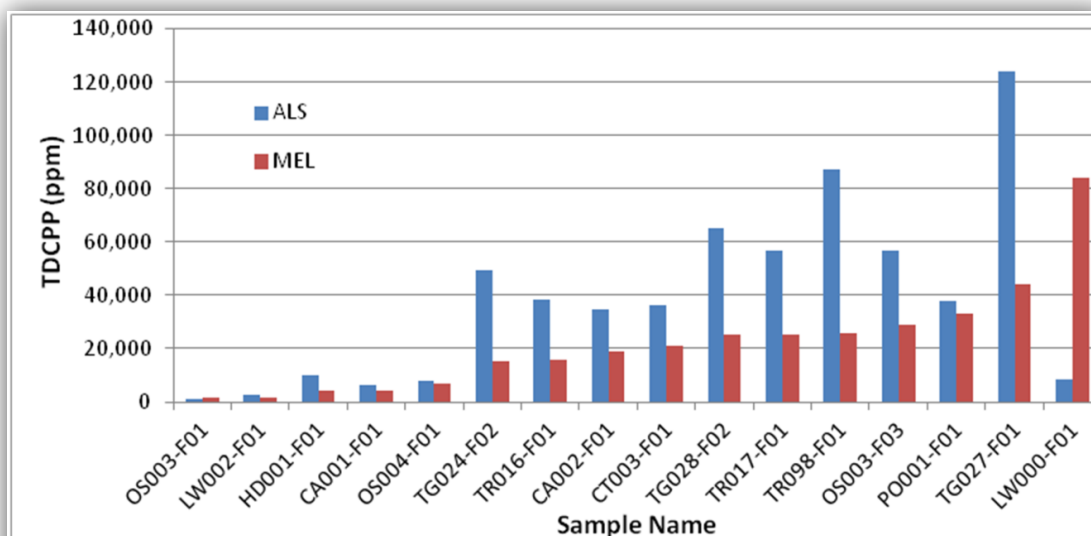


Figure 5: Comparison of laboratory levels in TDCPP in products.

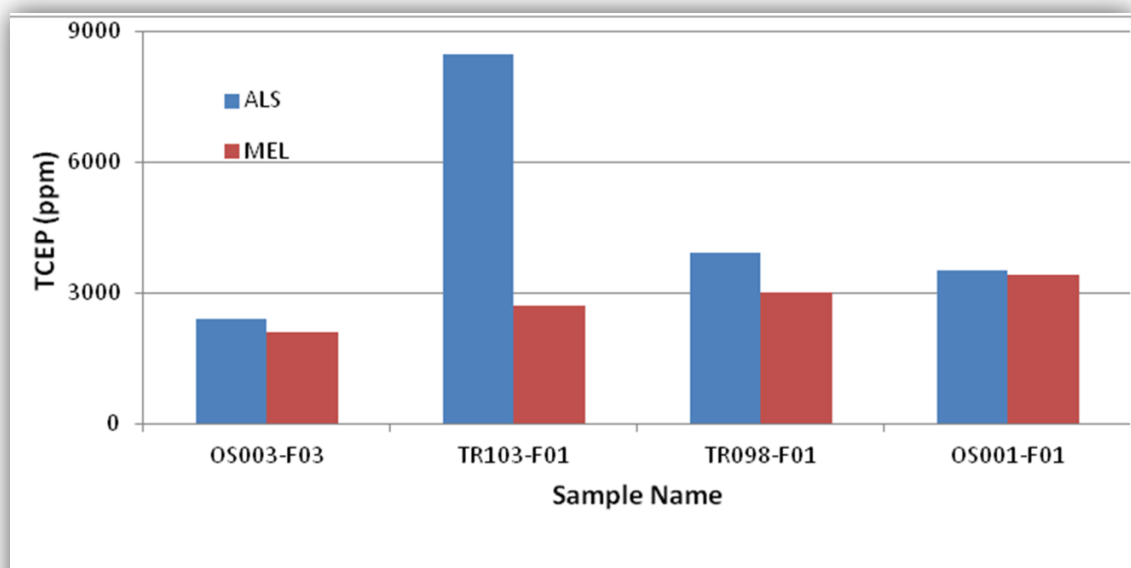


Figure 6: Comparison of laboratory levels of TCEP in products.

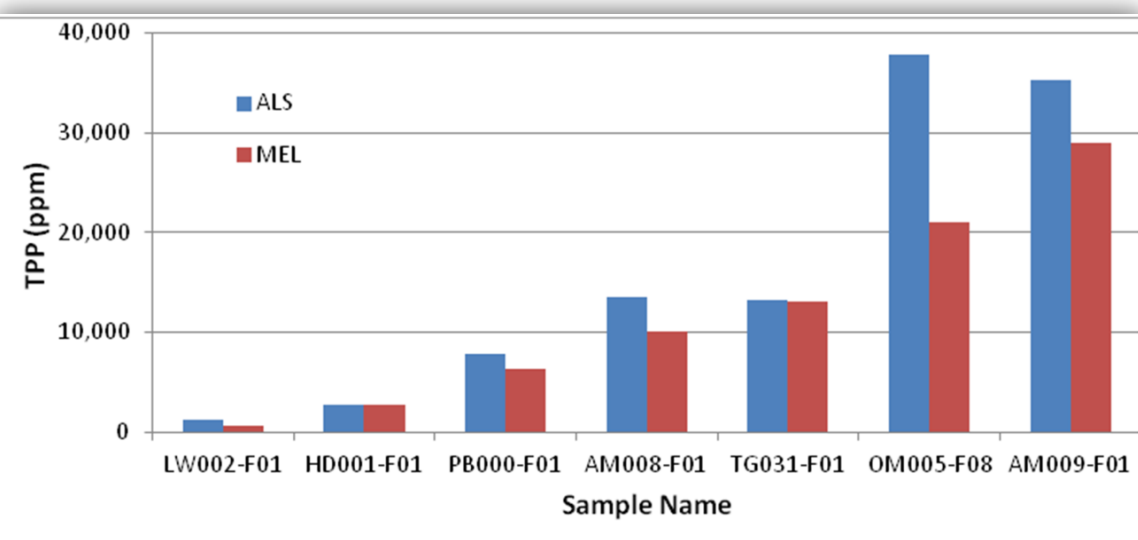


Figure 7: Comparison of laboratory levels of TPP in products.

Summary

Results from this study support the following conclusions:

- Most of the samples did not contain PBDEs above the 100 ppm reporting limit. Only three of 163 product components tested contained PBDEs above 1,000 ppm. Manufacturers have largely moved away from using PBDEs and their products are compliant with Washington regulations.
- Numerous products were found to contain chlorinated phosphate flame retardants and particularly TDCPP. A majority of these samples were foam and many were children's products. Several products contained multiple flame retardants suggesting mixtures are being used.
- A few samples contained TPP. Most were foam but some were fabric. XRF results indicate that three of these samples also contained bromine indicating a halogenated flame retardant was likely used in the mixture. These samples were children's products.
- Products were found to contain HBCD and TBBPA. These flame retardants were mainly found in plastic samples. Until companies reported as required under the CSPA, Ecology did not have information to indicate these flame retardants were used in children's products.
- V6 is being used in several children's products. Most of the samples were foam but one was fabric. TCEP is an impurity in V6 and levels in these products were above the CSPA criteria for TCEP.
- TCEP was found in seven samples from six products at levels above the CSPA reporting criteria. HBCD was found in one product at levels above the CSPA reporting criteria.
- While foam was the matrix where most of the analytes were detected, plastic had the most frequent detects via XRF for bromine. This suggests flame retardants not assessed in this study are in use.
- Problems exist with the analysis of both deca-PBD-ethane and RDP because of the inability to obtain and create adequate standards.
- Aside from deca-PBDethane and RDP, flame retardants can be measured in products at levels necessary to determine compliance with Washington regulations.

Recommendations

Based on findings of this study, the following recommendations are made:

- While not as effective as bromine, future sampling efforts for flame retardants should include a screening of chlorine using XRF. Reference standards should be purchased to ensure that the instrument is working within acceptable limits.
- Follow-up sampling and analysis should consider alternative analytes such as Firemaster® 550, commonly found in foam.
- Additional alternative brominated flame retardants should be assessed in plastics, as many plastic samples were found from screening analyses to contain bromine. For the majority of samples, however, the brominated alternatives were not identified.
- In order to assess TPP use, additional analytes commonly found in mixtures with TPP, such as Firemaster® 550 and tris(4-(tert-butyl)phenyl) phosphate should be measured.
- Children's products that meet the flammability requirements of CPAI-84 should be tested to understand what, if any, flame retardants are being used.
- Additional analysis of children's foam furniture should be performed. This group was most likely to have levels above the CSPA reporting criteria.
- Ecology should work with manufacturers and other interested parties to obtain pure standards for analysis.

References

- European Union. Risk Assessment Report of 2,2-bis(chloromethyl) trimethylene bis[bis(2-chloroethyl)phosphate] (V6). 2007. CAS No: 38051-10-4.
http://ec.europa.eu/health/ph_risk/committees/04_scher/docs/scher_o_067.pdf, accessed 11/2013.
- Allen, Joseph G., Michael D. McClean, Heather M. Stapleton and Thomas F. Webster, 2008. *Linking PBDEs in House Dust to Consumer Products using X-ray Fluorescence*, **Environ. Sci. Technol.**, 42, pages 4222-4228.
- California Environmental Protection Agency (Cal EPA), 2011. A chemical listed effective October 28, 2011 as known to the state of California to cause cancer TRIS (1,3-dichloro-2-propyl) Phosphate (TDCPP) (Cas No. 13674-87-8), available at: http://oehha.ca.gov/prop65/prop65_list/102811list.html, accessed 4/2012.
- Washington Department of Health (DOH), 2010. *Rationale for Reporting List of Chemicals of High Concern to Children Prepared by the Washington State Department of Health for the Children's Safe Product Act-Deca-BDE*, available at: www.ecy.wa.gov/programs/swfa/cspa/pdf/1163195.pdf, accessed 04/2012.
- Washington Department of Ecology (Ecology), 2006, Washington State Polybrominated Diphenyl Ether (PBDE) Chemical Action Plan: Final Plan, Ecology Publication No. 05-07-048, available at: <https://fortress.wa.gov/ecy/publications/SummaryPages/0507048.html>, accessed 3/2014.
- Ecology, 2008. Guidance Regarding Implementation of Chapter 70.76 RCW
www.ecy.wa.gov/programs/swfa/pbt/docs/DeMinimusPBDEGuidance.pdf
- Ecology, 2009. Alternatives to Deca-BDE in Televisions and Computers and Residential Upholstered Furniture, Publication No. 09-07-041, available at:
<https://fortress.wa.gov/ecy/publications/SummaryPages/0907041.html>, accessed 4/2012.
- Ecology, 2011a. The Reporting List of Chemicals of High Concern to Children (CHCC) website, available at: www.ecy.wa.gov/programs/swfa/cspa/chcc.html, accessed 3/2014.
- Ecology, 2011b. Chapter 173-334 WAC, Children's Safe Product Reporting Rule, available at: <https://fortress.wa.gov/ecy/publications/SummaryPages/173334.html>, accessed 4/2012.
- Ecology, 2012. Quality Assurance Project Plan: Flame Retardants in General Consumer and Children's Products, Publication number 12-07- 025, 29 pages,
<https://fortress.wa.gov/ecy/publications/SummaryPages/1207025.html>.
- Ecology, 2013. Addendum #1 to Quality Assurance Project Plan: Flame Retardants in General Consumer and Children's Products, Publication number 12-07-025A, 8 pages,
<https://fortress.wa.gov/ecy/publications/SummaryPages/1207025a.html>.
- Ecology, 2013. Children's Safe Product Act – Reporting Rule – WAC 173-334 Reporting Guidance – Practical Quantification Limits (PQLs)
www.ecy.wa.gov/programs/swfa/cspa/pdf/cspaguide_pql.pdf, accessed 3/2014.

Ecology, 2014. Product Sampling Standard Operating Procedure, Publication number 14-04-013, 11 pages, <https://fortress.wa.gov/ecy/publications/SummaryPages/1404013.html>.

EPA, 2010. DecaBDE Phase-out Initiative web page, available at: www.epa.gov/oppt/existingchemicals/pubs/actionplans/deccadbe.html, accessed 4/2012.

European Union (EU), 2009. Risk Assessment Report for Tris (2-chloroethyl) phosphate, TCEP, available at: http://esis.jrc.ec.europa.eu/doc/risk_assessment/REPORT/tcepreport068.pdf, accessed 4/2012.

Fang, M., T. F. Webster, D. Gooden, E. Cooper, M. McClean, C. Carignan, C. Makey & H. Stapleton, 2013. *Investigating a Novel Flame Retardant Known as V6: Measurements in Baby Products, House Dust, and Car Dust*, **Environ. Sci. Technol.**, 47, pages 4449-4454.

Keller, A.S., P.R. Nikhilesh, T.F. Webster, H.M. Stapleton. 2014. *Flame retardant applications in camping tents and potential exposure*. Environmental Science & Technology Letters. <http://pubs.acs.org/doi/abs/10.1021/ez400185y> last accessed 3/13/14.

Kim, J. W., T. Isobe, K. H. Chang, A. Amano, R. H. Maneja, P. B. Zamora, F. P. Siringan, & S. Tanabe, 2011. *Levels and distribution of organophosphorus flame retardants and plasticizers in fishes from Manila Bay, the Philippines*, **Environ Pollut** 159, pages 3653-3659.

Kuehl D. W. & R. Haebler, 1995. *Organochlorine, organobromine, metal, and selenium residues in bottlenose dolphins (Tursiops truncatus) collected during an unusual mortality event in the Gulf of Mexico, 1990*, **Arch Environ Contam Toxicol**, 28, pages 494-499.

Loraine G. A. & M. E. Pettigrove, 2006. *Seasonal variations in concentrations of pharmaceuticals and personal care products in drinking water and reclaimed wastewater in southern California*. **Environ Sci Technol**, 40, pages 687-695.

Marklund, A., 2005. Levels and Sources of Organophosphorus Flame Retardants and Plasticizers in Indoor and Outdoor Environments, Umea University, Department of Chemistry, 69 pages.

National Industrial Chemicals Notification and Assessment Scheme (NICNAS), 2001. Triphosphates Priority Existing Chemical Assessment Report No. 17, Commonwealth of Australia, available at: www.nicnas.gov.au/Publications/CAR/PEC/PEC17/PEC_17_Full_Report_PDF.pdf, accessed 4/2012.

Organization for Economic and Cooperative Development (OECD), 2012 (actual date of publication not provided). Safety Information Data Sheet (SIDS) for Tris (1-chloro-2-phosphphate) phosphate, CAS No. 16374-84-5, available at: www.chem.unep.ch/irptc/sids/oecdsids/13674845.pdf, accessed 4/2012.

Ross, Peter S., 2006. *Fireproof killer whales (orcinus orca): flame-retardant chemicals and the conservation imperative in the charismatic icon of British Columbia, Canada*, **Can. J. Fish Aquat. Sci.**, 63, pages 224-234.

Stapleton, Heather M. Joseph G. Allen, Shannon M. Kelly, Alex Konstantinov, Susan Klosterhaus, Deborah Watkins, Michael D. McClean and Thomas F. Webster, 2008. *Alternate and New Brominated Flame Retardants Detected in U.S. House Dust*, **Environ. Sci. Technol.**, 42, pages 6910-6916.

Stapleton, Heather M., Susan Klosterhaus, Sarah Eagle, Jennifer Fuh, John D. Meeker, Arelene Blum and Thomas F. Webster, 2009. *Detection of Organophosphate Flame Retardants in Furniture Foam and U.S. House Dust*, **Environ. Sci. Technol.**, 43, pages 7490-7495.

Stapleton, Heather, Susan Klosterhaus, Alex Keller, P. Lee Ferguson, Saskia van Bergen, Ellen Cooper, Thomas F. Webster and Arlene Blum, 2011. *Identification of Flame Retardants in Polurethane Foam Collected from Baby Products*, **Environ. Sci. Technol.**, 45, pages 5323-5331.

Stapleton, Heather, Smriti Sharma, Gordon Getzinger, P. Lee Ferguson, Michelle Gabriel, Thomas F. Webster and Arlene Blum, 2012. *Novel and High Volume Use Flame Retardants in US Couches Reflective of the 2005 PentaBDE Phase Out*, **Environ. Sci. Technol.**, 46, pages 13432-13439.

Stone, Alex and Damon Delistraty, 2010, 'Sources of toxicity and exposure information for identifying chemicals of high concern to children', **Env. Impact Assess. Review**, 30, pages 380–387, available at: www.sciencedirect.com/science/article/pii/S0195925509001437.

Sundkvist, Anneli Marklund, Ulrika Olofsson and Peter Haglund, 2010. *Organophosphorus flame retardants and plasticizers in marine and fesh water biota and in human milk*, **J. Environ. Monit.**, 12, pages 943-951.

Van den Eede, N., A. C. Dirtu, H. Neels, A. Covaci, 2011. *Analytical developments and preliminary assessment of human exposure to organophosphate flame retardants from indoor dust*. **Environ Int** 37, pages 454-461.

Washington Department of Health (DOH), 2010. *Rationale for Reporting List of Chemicals of High Concern to Children Prepared by the Washington State Department of Health for the Children's Safe Product Act-Deca-BDE*, available at: www.ecy.wa.gov/programs/swfa/cspa/pdf/1163195.pdf, accessed 04/2012.

Washington State Legislature (Washington), 2007. Chapter 70.76 RCW Polybrominated diphenyl ethers-flame retardants, available at: <http://apps.leg.wa.gov/RCW/default.aspx?cite=70.76&full=true>, accessed 4/2012.

Washington, 2008. Chapter 70.240 Children's Safe Products, available at: <http://apps.leg.wa.gov/RCW/default.aspx?cite=70.240&full=true>, accessed 4

Acronyms, Abbreviations, and Appendices

Acronyms

CHCC	Chemical of high concern to children
CSPA	Washington's Children's Safe Product Act
DOH	Washington State Department of Health
Ecology	Washington State Department of Ecology
EPA	U.S. Environmental Protection Agency
FDA	Food and Drug Administration
HBCD	Hexabromocyclododecane
FM550	Firemaster® 550
LCS	Laboratory control sample
MEL	Manchester Environmental Laboratory
MQO	Measurement quality objective
OECD	Organization of Economic and Cooperative Development
RCW	Revised Code of Washington
RDP	Resorcinol diphenyl phosphate
PBDEs	Polybrominated diphenyl ethers
PBD-ethanes	Polybrominated diphenyl ethanes
TBBPA	Tetrabromobisphenol-A
TCEP	Tris(2-chloroethyl) phosphate
TCPP	Tris(1-chloro-2-propyl) phosphate
TDCPP	Tris(1,3-dichloro-2-propyl) phosphate
TPP	Triphenyl phosphate
U-OPFR	Unidentified chlorinated organophosphate flame retardant
V6	2,2-bis(chloromethyl)-propane-1,3-diyltetrakis(2-chloroethyl) bisphosphate
WAC	Washington Administrative Code
XRF	X-ray fluorescence

Data Flags

- J The analyte was positively identified; the associated numerical value is the approximate concentration of the analyte in the sample.
- UJ The analyte was analyzed for and was not present above the level of the associated value. The associated numerical value may not accurately or precisely represent the concentration necessary to detect the analyte in this sample.
- U The analyte was not detected at or above the reported sample quantitation limit.

Units of Measurement

ppm part per million

Appendix 1. Information on the components sent to the Laboratory.

<i>Bolded green components were sent to both Laboratories.</i>			<i>Blue highlighted results are children's products.</i>			
Component ID	Material	Segment	Product Description	Component	Flame retardant detected in order of relative abundance	XRF Br > 1%
SK023-F03	Fabric	Toys/Games	Baby Activity Gym	orange fabric covering all of base	-	~.5%
SK023-F06	Fabric	Toys/Games	Baby Activity Gym	animal fabric on top of base of gym	-	
SK025-F01	Fabric	Toys/Games	Baby Activity Gym	red fabric on bottom of mat	-	
SK025-F03	Fabric	Toys/Games	Baby Activity Gym	patterned fabric on top of mat	-	
SK024-F03	Plastic	Baby Care	Baby Carrier	black mesh fabric	-	
TR103-F01	Fabric	Baby Care	Baby Carrier	light grey fabric	V6, TCEP, TCPP, TDCPP	
TR103-F12	Fabric	Baby Care	Baby Carrier	composite of light and dark grey fabric	TCEP, TCPP, TDCPP	
TR019-F02	Stuffing	Baby Care	Baby neck support	stuffing	-	
TR020-F02	Stuffing	Baby Care	Baby neck support	stuffing	-	
WM095-F01	Styrofoam	Household/Office Furniture/Furnishings	Bean Bag Chair	poly foam beans	HBCD	+
TR018-F01	Foam	Household/Office Furniture/Furnishings	Belly Pillow	inner foam	-	
TG027-F01	Foam	Baby Care	Booster Seat	interior foam	TDCPP, TCPP, TCEP	
TG033-F01	Stuffing	Household/Office Furniture/Furnishings	Booster Seat Pad	padding	-	
TR105-F06	Fabric	Baby Care	Car Seat	polka dotted fabric in seat	-	~.5%
TR105-F09	Fabric	Baby Care	Car Seat	composite of green and brown fabric from head support	-	+
TR105-F11	Fabric	Baby Care	Car Seat	dark brown fabric in seat	-	+
WM097-F01	Foam	Baby Care	Car Seat Handle Cushion	gray exterior foamy fabric	-	
DT042-F02	Foam	Baby Care	Changing mat	foam	-	
TG024-F02	Foam	Baby Care	Changing mat	foam	TDCPP, TCPP, TCEP	
TG028-F02	Foam	Baby Care	Changing Pad	interior foam	TDCPP, TCPP, TCEP	
TR015-F01	Foam	Baby Care	Changing Pad	interior foam	TCPP, TDCPP, TCEP	
MR000-F03	Foam	Household/Office Furniture/Furnishings	Mattress	top foam under batting	-	
MR000-F06	Fabric	Household/Office Furniture/Furnishings	Mattress	felt fabric under bottom fabric	-	
CT006-F01	Fabric	Baby Care	Mattress	cover fabric shiny side down	-	

<i>Bolded green components were sent to both Laboratories.</i>			<i>Blue highlighted results are children's products.</i>			
Component ID	Material	Segment	Product Description	Component	Flame retardant detected in order of relative abundance	XRF Br > 1%
CT006-F03	Stuffing	Baby Care	Mattress	internal batting	-	
CT006-F04	Foam	Baby Care	Mattress	top foam layer	-	
CT006-F05	Foam	Baby Care	Mattress	internal thick foam	-	
TG023-F01	Fabric	Clothing	Pajamas	grey pjs	-	
TG051-F01	Fabric	Clothing	Pajamas	PJ Fabric	-	
WM065-F01	Fabric	Clothing	Pajamas	PJ Pant Fabric	-	
CP000-F01	Fabric	Clothing	Pajamas	Fleece footie pjs	-	
FM040-F01	Fabric	Clothing	Bathrobe	Bathrobe Fabric	-	
WM063-F01	Fabric	Clothing	Bathrobe	Bathrobe fabric	-	
AM007-F01	Foam	Household/Office Furniture/Furnishings	Chair	arm rest foam	-	
AM008-F01	Foam	Household/Office Furniture/Furnishings	Chair	Foam	TPP	+
AM009-F01	Foam	Household/Office Furniture/Furnishings	Chair	Foam	TPP	
OS000-F02	Foam	Household/Office Furniture/Furnishings	Chair	blue/green seat cushion foam	-	
OS000-F04	Foam	Household/Office Furniture/Furnishings	Chair	pink body foam	-	
OS001-F01	Foam	Household/Office Furniture/Furnishings	Chair	foam	V6, TCPP, TCEP, TDCPP	
OS002-F01	Foam	Household/Office Furniture/Furnishings	Chair	upper portion inner foam	TDCPP, TCPP, TCEP	
OS003-F01	Foam	Household/Office Furniture/Furnishings	Chair	head cushion fluff	TDCPP, TCPP, TCEP	
OS003-F03	Foam	Household/Office Furniture/Furnishings	Chair	head cushion fluff yellow foam chunks	TDCPP, TCPP, TCEP	
PB000-F01	Foam	Household/Office Furniture/Furnishings	Chair	foam	TPP	+
TG031-F01	Foam	Household/Office Furniture/Furnishings	Chair	interior foam	TPP	+
TR097-F01	Foam	Household/Office Furniture/Furnishings	Chair	chair foam	-	
TR098-F01	Foam	Household/Office Furniture/Furnishings	Chair	chair foam	TDCPP, V6, TCPP, TCEP	
WM094-F01	Foam	Household/Office Furniture/Furnishings	Chair	cushion foam	TCPP, TDCPP, TCEP	
OS004-F01	Foam	Household/Office Furniture/Furnishings	Chair Accessory	foam from cushion	TDCPP, TCPP, TCEP	
CT000-F01	Stuffing	Household/Office Furniture/Furnishings	Pillow	stuffing	-	

<i>Bolded green components were sent to both Laboratories.</i>			<i>Blue highlighted results are children's products.</i>			
Component ID	Material	Segment	Product Description	Component	Flame retardant detected in order of relative abundance	XRF Br > 1%
TG029-F02	Stuffing	Household/Office Furniture/Furnishings	Pillow	interior stuffing	-	
TG030-F01	Stuffing	Household/Office Furniture/Furnishings	Pillow	inner stuffing	-	
BY002-F01	Plastic	Toys/Games	Tablet	Tablet- for kid- plug		+
BY002-F08	Plastic	Toys/Games	Tablet	Tablet- for kid- green plastic around battery	PBD-ethane	+
BY002-F11	Plastic	Toys/Games	Tablet	Tablet- for kid- motherboard	-	~.5%
TR017-F01	Foam	Baby Care	crib wedge	interior foam	TDCPP, TCPP, TCEP	
TR021-F02	Foam	Baby Care	Pee protector	interior foam	TDCPP, TCPP, TCEP	
TR022-F02	Foam	Baby Care	Pee protector	Foam	-	
CT001-F01	Stuffing	Toys/Games	Play Pad	Stuffing	-	
TR016-F01	Foam	Baby Care	Portable Crib Pad	interior foam	TDCPP, TCPP, TCEP	
TG026-F03	Stuffing	Toys/Games	Security Blanket	owl stuffing	-	
TG025-F02	Stuffing	Toys/Games	Stuffed Animal	stuffing	-	
DT037-F02	Stuffing	Toys/Games	Stuffed Toy	stuffing	-	
DT039-F01	Stuffing	Toys/Games	Stuffed Animal	Interior stuffing	-	
FM042-F01	Plastic	Beauty/Personal Care/Hygiene	Hair Dryer	Side of dryer	-	
BL008-F01	Plastic	Electrical Supplies	6 Outlet surge protector	Body of strip	-	
BL036-F04	Plastic	Computing	7 inch internet tablet	7 inch internet tablet- motherboard	-	+
CT002-F01	Foam	Household/Office Furniture/Furnishings	Bath Mat	Foam	-	
FM038-F01	Plastic	Electrical Supplies	Battery Charger	Charger body	-	~.5%
OM000-F02	Plastic	Electrical Supplies	Battery Charger	Silver body	TBBPA, TPP	+
WM066-F01	Plastic	Electrical Supplies	Battery Charger	Charger body	TPP	
CT010-F02	Plastic	Audio Visual/Photography	blue ray disc/dvd player	plastic exterior	-	
CT010-F05	Plastic	Audio Visual/Photography	blue ray disc/dvd player	white plastic	TPP	
CT010-F08	Plastic	Audio Visual/Photography	blue ray disc/dvd player	larger motherboard	-	+
WS002-F01	Plastic	Automotive	Car charger	Side	TBBPA	+

<i>Bolded green components were sent to both Laboratories.</i>			<i>Blue highlighted results are children's products.</i>			
Component ID	Material	Segment	Product Description	Component	Flame retardant detected in order of relative abundance	XRF Br > 1%
HD000-F01	Foam	Building Products	Carpet padding	Carpet padding	TDCPP, TPP, TCP, TCEP	
HD001-F01	Foam	Building Products	Carpet padding	Carpet padding with red backing	TDCPP, TPP, TCP, TCEP	
LW000-F01	Foam	Building Products	Carpet padding	Carpet padding	TDCPP, TPP, TCP, TCEP, V6	
LW001-F01	Foam	Building Products	Carpet padding	Carpet padding	TDCPP, TCP, TPP, TCEP	
LW002-F01	Foam	Building Products	Carpet padding	Carpet padding	TDCPP, TCP, TPP, V6, TCEP	
TG056-F01	Plastic	Home Appliances	Ceramic Heater	Space heater side plastic	-	+
CT003-F01	Foam	Household/Office Furniture/Furnishings	Chair Pad	Foam	TDCPP, TCP, TPP, TCEP	
TG050-F01	Plastic	Communications	Charger	Plug in	-	+
FM043-F01	Plastic	Home Appliances	Clothing Iron	Silver plastic closest to heating element	TBBPA	+
FM043-F02	Plastic	Home Appliances	Clothing Iron	Cord	-	
WM069-F01	Plastic	Home Appliances	Coffee Pot	Bottom plastic closest to heating element	-	+
FM039-F01	Plastic	Electrical Supplies	compact power station	Charger body	-	
BB000-F02	Plastic	Home Appliances	Dehydrator	grey plastic around top	-	
BB000-F03	Plastic	Home Appliances	Dehydrator	grey plastic around top	-	
BB000-F04	Plastic	Home Appliances	Dehydrator	grey plastic above coil	-	+
GR008-F01	Plastic	Electrical Supplies	Electrical Supply		-	
FM041-F01	Plastic	Home Appliances	Fan Forced Heater	Space heater side plastic	-	+
WM070-F01	Plastic	Home Appliances	Fan-Forced Heater	Space heater side plastic	-	
GR014-F01	Oil/Solvent	Home/Business Safety/Security/Surveillance	Fire Retardant spray	liquid spray	-	
GR015-F01	Plastic	Lawn/Garden Supplies	Flame resistant drainage tarp	Flame resistant drainage tarp plastic	-	+
GR010-F01	Plastic	clothing	Flame resistant lab coat	blue lab coat fabric	-	
GR013-F02	Plastic	Electrical Supplies	flame resistant sleeve	blue elastic fabric	-	
GR012-F01	Plastic	Home/Business Safety/Security/Surveillance	Flame retardant pad	blue pad fabric	-	
GR012-F02	Stuffing	Home/Business Safety/Security/Surveillance	Flame retardant pad	Stuffing	-	
TG053-F01	Plastic	Communications	Flip Phone	Phone-battery side	-	

<i>Bolded green components were sent to both Laboratories.</i>			<i>Blue highlighted results are children's products.</i>			
Component ID	Material	Segment	Product Description	Component	Flame retardant detected in order of relative abundance	XRF Br > 1%
HD002-F01	Foam	Camping	foam camping pad	foam camping pad	-	
CT007-F01	Fabric	Household/Office Furniture/Furnishings	Foam Pillow	external cloth cover	-	
CT007-F02	Fabric	Household/Office Furniture/Furnishings	Foam Pillow	internal cloth protector	-	
CT007-F03	Foam	Household/Office Furniture/Furnishings	Foam Pillow	internal gel blue beads	-	
CT007-F04	Foam	Household/Office Furniture/Furnishings	Foam Pillow	internal white foam	-	
TG052-F01	Plastic	Beauty/Personal Care/Hygiene	Hair Dryer	Side of dryer	-	
WM067-F01	Plastic	Beauty/Personal Care/Hygiene	Hair Dryer	Purple outer plastic	-	
WM067-F03	Plastic	Beauty/Personal Care/Hygiene	Hair Dryer	cord	-	
BL005-F01	Plastic	Beauty/Personal Care/Hygiene	Hair Flat Iron	Outside of flat iron	TBBPA	+
BL005-F03	Plastic	Beauty/Personal Care/Hygiene	Hair Flat Iron	Inner plastic - near element	-	
BL003-F01	Plastic	Home Appliances	Hand Mixer	Side of mixer	-	
WM068-F01	Plastic	Household/Office Furniture/Furnishings	Heated mattress foot warming pad	Heat pad	-	
WM068-F02	Plastic	Household/Office Furniture/Furnishings	Heated mattress foot warming pad	Top of controller	-	+
WM068-F03	Plastic	Household/Office Furniture/Furnishings	Heated mattress foot warming pad	Bottom of controller	TBBPA	+
WM068-F04	Plastic	Household/Office Furniture/Furnishings	Heated mattress foot warming pad	Connector	-	+
BL009-F01	Fabric	Household/Office Furniture/Furnishings	Heated Throw	Blanket Fabric	-	
BL009-F03	Plastic	Household/Office Furniture/Furnishings	Heated Throw	Large cord	-	
HD004-F01	Fiberglas	Building Products	Insulation	Fiberglas insulation	-	
BB001-F01	Fabric	Home Appliances	Ironing Board Cover	beige cover	-	
GR006-F01	Fabric	Home Appliances	Ironing Board Cover	exterior fabric	-	
OM005-F08	Plastic	Computing	LCD Monitor	largest motherboard - green one side tan the other	TPP	+
OM005-F09	Plastic	Computing	LCD Monitor	medium motherboard - green both sides	-	+
OM005-F10	Plastic	Computing	LCD Monitor	plastic back	-	
CT009-F03	Plastic	Audio Visual/Photography	LED TV	plastic covering power input	-	+
CT009-F08	Plastic	Audio Visual/Photography	LED TV	power supply	-	

Bolded green components were sent to both Laboratories.

Blue highlighted results are children's products.

Component ID	Material	Segment	Product Description	Component	Flame retardant detected in order of relative abundance	XRF Br > 1%
CT009-F09	Plastic	Audio Visual/Photography	LED TV	cable connector	-	
CT009-F12	Plastic	Audio Visual/Photography	LED TV	Mother board from TV	-	+
BL037-F02	Foam	Household/Office Furniture/Furnishings	Mattress	top batting	-	
CT008-F01	Foam	Household/Office Furniture	Mattress Pad	Foam	-	
DT044-F02	Fabric	Clothing	Oven Mitt		-	
GR005-F01	Fabric	Clothing	Oven Mitt	oven mitt	-	
BL004-F01	Plastic	Home Appliances	Personal Heater Fan	Space heater side plastic	-	+
OR000-F01	Plastic	Storage/Haulage Containers	Plastic Pallet	structural foam pallet black	-	
BC000-F01	Plastic	Storage/Haulage Containers	Plastic Pallet	Plastic	-	
GR000-F01	Plastic	Storage/Haulage Containers	Plastic Pallet	Black	-	
GR001-F01	Plastic	Storage/Haulage Containers	Plastic Pallet	recycled plastic pallet	-	
UL000-F01	Plastic	Storage/Haulage Containers	Plastic Pallet	Plastic Pallet	-	
WS001-F01	Plastic	Camping	Pot and Pan Support	Camp stove converter-Orange Plastic stand	-	
HD005-F01	Plastic	Home Appliances	Power Cord	cord casing	-	
BY000-F01	Plastic	Computing	Projector	Projector-Top plastic casing	-	
WS006-F01	Plastic	Camping	Propane Heater	Side Casing	-	+
GR002-F01	Fabric	Clothing	Protective Clothing	yellow fabric	-	
GR003-F01	Plastic	Clothing	Protective Glove	glove rubber palm fabric	HBCD	
WW000-F01	Plastic	Clothing	Protective Gloves	black rubber surface	-	
WW001-F01	Plastic	Clothing	Protective Gloves	Vinyl glove	-	
WS004-F01	Plastic	Automotive	Quick Jumper Cables	Yellow Casing	-	
LW004-F01	Plastic	Home Appliances	range receptacle	big black plastic connector	-	+
DS000-F01	Fabric	Clothing	Reflective Clothing	Orange Fabric	-	
DS000-F02	Plastic	Clothing	Reflective Clothing		-	
MW000-F01	Fabric	Clothing	Reflective Vest	Orange Fabric	-	
MW000-F02	Fabric	Clothing	Reflective Vest	yellow fabric	-	

<i>Bolded green components were sent to both Laboratories.</i>			<i>Blue highlighted results are children's products.</i>			
Component ID	Material	Segment	Product Description	Component	Flame retardant detected in order of relative abundance	XRF Br > 1%
SG000-F01	Plastic	Clothing	Reflective Vest	Reflective tape	-	
SG000-F02	Fabric	Clothing	Reflective Vest	Fabric	-	
PO001-F01	Foam	Household/Office Furniture/Furnishings	Seat Cushion	red fabric	TDCPP, TCPP, TCEP	
OM003-F01	Plastic	Stationery/Office Machinery/Occasion Supplies	Shredder	Large grey plastic holder piece	TPP	+
OM003-F02	Plastic	Stationery/Office Machinery/Occasion Supplies	Shredder	top plastic	TBBPA, TPP	+
CA002-F01	Fabric	Camping	Tent	Tent- black tent fabric	TDCPP, V6, TCEP	
CA002-F02	Fabric	Camping	Tent	Tent- silver tent fabric	TDCPP, TCEP, TCPP	
CA001-F01	Fabric	Camping	Tent floor liner	Tent floor liner- green fabric	TDCPP, TCEP	
FM044-F01	Plastic	Home Appliances	Toaster	Plastic ends of toaster	-	
TG057-F01	Plastic	Home Appliances	Toaster	Side of plastic toaster	-	
WM071-F01	Plastic	Home Appliances	Toaster	Side of plastic toaster	-	
CA000-F01	Fabric	Camping	Tube Tent	Tube Tent- Emergency Shelter- orange plastic	-	
RS001-F01	Plastic	Computing	Wireless Router	Wireless Router-Case	-	

Appendix 2

Components containing target analytes above 1000 ppm from ALS or MEL

Component ID	Analyte	ALS result (ppm)	MEL result (ppm)	% RPD	Material
CA001-F01	TDCPP	6190	4300	36.0	Fabric
CA002-F01	TDCPP	34900	19000 J	59.0	Fabric
CT003-F01	TDCPP	36200	21000 J	53.1	Foam
HD001-F01	TDCPP	9880 J	4100 J	82.7	Foam
LW000-F01	TDCPP	8430	84000 J	163.5	Foam
LW002-F01	TDCPP	2550	1800 J	34.5	Foam
OS003-F03	TDCPP	56600	29000 J	64.5	Foam
OS004-F01	TDCPP	8000	7000 J	13.3	Foam
PO001-F01	TDCPP	37900	33000	13.8	Foam
TG024-F02	TDCPP	49500	15000 J	107.0	Foam
TG027-F01	TDCPP	124000	44000 J	95.2	Foam
TG028-F02	TDCPP	65400	25000 J	89.4	Foam
TR016-F01	TDCPP	38200 J	16000 J	81.9	Foam
TR017-F01	TDCPP	56800	25000 J	77.8	Foam
TR098-F01	TDCPP	87100	26000	108.0	Foam
OS003-F01	TDCPP	926 J	1600	53.4	Foam
OS001-F01	TCEP	3530	3400 J	3.8	Foam
OS003-F03	TCEP	2390	2100 J	12.9	Foam
TR098-F01	TCEP	3930	3000 J	26.8	Foam
TR103-F01	TCEP	8490	2700	103.5	Foam
CT003-F01	TCPP	3720	3000	21.4	Foam
LW000-F01	TCPP	1040	690	40.5	Foam
OS001-F01	TCPP	20900	12000	54.1	Foam
OS003-F03	TCPP	16600	12000 J	32.2	Foam
PO001-F01	TCPP	6830	7000 J	2.5	Foam
TG024-F02	TCPP	3200	2800 J	13.3	Foam
TG027-F01	TCPP	43200	15000	96.9	Foam
TG028-F02	TCPP	3990	4900 J	20.5	Foam
TR016-F01	TCPP	29300	9700	100.5	Foam
TR017-F01	TCPP	8230	7700	6.7	Foam
TR098-F01	TCPP	19300	16000	18.7	Foam
TR103-F01	TCPP	1690	640 J	90.1	Foam
AM008-F01	TPP	13500	10000 J	29.8	Foam
AM009-F01	TPP	35300	29000 J	19.6	Foam
HD001-F01	TPP	2710 J	2700	0.4	Foam
LW002-F01	TPP	1150	560	69.0	Foam
OM005-F08	TPP	37800	21000	57.1	Plastic
PB000-F01	TPP	7770	6300	20.9	Foam
TG031-F01	TPP	13200	13000 J	1.5	Foam